

PURSUIT[®] XL

ENGINEERED FOR POWERFUL PERFORMANCE AND EXCEPTIONAL STYLE

PURSUIT[®] XL

The Pursuit[®] XL is an exciting blend of power and precision that features full suspension, low-profile tires and a hydraulic-sealed brake system. Standard luxury touches like a wraparound delta tiller and LED lighting make it the ultimate outdoor scooter.



Delta tiller for optimal handling



Pillow top seat with sliders



LED lighting package



Large 13" pneumatic tires




Full suspension



Pride[®]
Mobility Products Corp.
— LIVE YOUR BEST™ —

Pursuit® XL


Model number:	SC714
Weight capacity:	400 lbs.
Maximum speed ¹ :	Up to 8.25 mph high speed
Ground clearance ² :	5.25"
Turning radius:	82.5"
Overall length ² :	56.25"
Overall width ² :	27.25"
Front tires:	13" pneumatic
Rear tires:	13" pneumatic
Seat-to-deck:	27"-28.5"
Suspension:	Front & Rear
Range per charge ¹ :	28.94 miles with group 27 (100 AH) 21.70 miles with group 24 (70 AH)
Total weight without batteries:	205.5 lbs. (main frame)
Standard seating:	Type: Reclining, pillow top, high-back seat with headrest Weight: 48 lbs. Material: Black Vinyl Dimensions: width - 19.5" (overall width) depth - 17.5" (effective depth)
Drive system:	Rear-wheel drive, sealed transaxle, 24-volt DC motor
Dual braking system:	Regenerative and electromechanical
Battery requirements ³ :	Two 12-volt, deep cycle Size: Group 27 (100 AH) Group 24 (70 AH) Weight ⁴ : Group 24: 55 lbs. each; Group 27: 60.5 lbs.
Battery charger:	8-amp, off-board
Colors:	

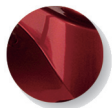

Candy Apple Red



Viper Blue

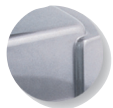

Silver

¹Varies with user weight, terrain type, battery amp-hour (AH), battery charge, battery condition and tire condition.

²Due to manufacturing tolerances and continual product improvement, this specification can be subject to a variance of + or - 3%.

³AGM or Gel-Cell type required.

⁴Battery weight may vary based on manufacturer.

⁵Includes standard seat.

NOTE: All specifications subject to change without notice.

Features

- Infinitely-adjustable tiller angle with conveniently located handle
- Powerful drivetrain; features 180 amp microprocessor-based controller
- Deluxe reclining high-back seat with headrest and sliders
- Full lighting package that includes directional signals and lower headlight
- Delta tiller with wraparound handles
- Large 13" low-profile, pneumatic tires for excellent outdoor performance
- Front and rear suspension
- Easily accessible tie down points (for transport of unoccupied scooter)
- Front basket standard
- Rearview mirror standard
- Hydraulic handbrake

Options

- Rear basket
- Safety flag
- Durable weather cover



The information contained herein is correct at the time of publication; we reserve the right to alter specifications without prior notice. Speed & range vary with user weight, terrain type, battery charge, battery condition and tire pressure. Due to manufacturing tolerances and continual product improvement, length, width, turning radius and ground clearance can be subject to variance of + or - 3%.

800-800-8586 (US) • 888-570-1113 (Canada) • www.pridemobility.com

