

Buzzaround 3-wheel GB118 4-wheel GB148



The Golden Buzzaround Extreme Scooter features the portability and easy disassembly of a lightweight scooter with the amazing benefits of a full-size scooter! Standard comfort features on the Buzzaround EX include a stylish guilted high back seat, an infinite adjustable tiller at your fingertips, and a USB charging port right under the dashboard. The 4-Wheel model features a sporty new LED V-shaped running light on the front for maximum safety and visibility. The Buzzaround EX also comes standard with 2-U1 batteries - industry leading and the largest battery in this category - for maximum operating range up to 18 miles. With a 350 lb. weight capacity, comfort spring suspension, and 40.5" turning radius, the Buzzaround EX offers a smooth ride. Go the extra mile on a BuzzAround Extreme Scooter!



Standard Colors:







Coastal Blue (Included in Box)

NEW Sporty Optional Colors:



Crisp White



Canary Yellow





Stone Grey



Army Green (Matte)



Desert Tan

Sporty New LED Running Light on the 4-W/ EX!

Standard Features:

- 350 lb. weight capacity
- 4" ground clearance
- Full size & portable
- Larger 90 amp controller
- Larger 320 watt motor
- Larger 9" tires
- U1 battery for maximum range
- Ultrabright LED headlight & tail light
- Comfort suspension springs
- Full size, yet portable & easy to disassemble




Model	GB118	GB148
Description	EX	EX
Medicare Code	Not Coded	Not Coded
Weight Capacity	350 lbs.	350 lbs.
Drive Wheels	Rear	Rear
Maximum Speed	5 mph	5 mph
Operating Range up to ¹	18 miles	18 miles
Under Deck Clearance	4"	4"
Turning Radius	40.5"	54"
Type Batteries	2-U1	2-U1
Colors	Red, Blue	Red, Blue
Freewheel Mode	Yes	Yes
Electronic Speed Control	Yes	Yes
Electro-Mechanical Brakes	Yes	Yes
Charger	Off-board	Off-board
•		
Controller	90 Amp Dynamic R-Series	90 Amp Dynamic R-Series
Dimensions:		
Length	45"	44"
Width	22"	22"
Height (ground to top of back of seat)	35"-37"	35"-37"
Ground to Top of Seat	22.75"-24.75"	22.75"-24.75"
	6.75"	6.75"
Ground to Top of Deck	16"-18"	
Deck to Top of Seat		16"-18"
Number of Seat Height Positi		3
Size of Increments	1"	1"
Front Axle to Rear Axle	34.7"	33.8"
Weight of Unit: (assembled)	154 lbs.	161 lbs.
Front Half	46 lbs.	53 lbs.
Rear Half	36 lbs.	36 lbs.
Seat w/Arms	23 lbs.	23 lbs.
Batteries	24 lbs. per pack	24 lbs. per pack
Tires:	2 1 1001 po. pacit	2 1 1501 po. pas.
Standard Tires	Solid, low profile	Solid, low profile
Optional Tires	N/A	N/A
Front	9"	9"
Rear	9"	9"
	2"	2"
Rear Anti-Tip Wheels Standard Seat	Highback	– – – Highback
Seat Back Height	· ·	
- no headrest	18"	18"
Width x Depth	18" x 16"	18" x 16"
Color	Black Vinyl	Black Vinyl
Optional Seat:	Stadium Style	Stadium Style
Back Height - no headrest	16"	16"
Width x Depth	20"x 18"	20"x 18"
Color		
	Black Vinyl	Black Vinyl
Warranty:	1 :f_+: \ \\/	1 :f_+: \ \\/
Frame	Lifetime Warranty	Lifetime Warranty
D: T:		
Drive Train	2 Year Warranty	
Drive Train Electronics ² Batteries	2 Year Warranty 2 Year Warranty 1 Year Warranty	2 Year Warranty 2 Year Warranty 1 Year Warranty

^{1.} Battery range at 200 lbs. and will vary due to rider weight, drive surface, terrain, and battery type.

FDA Class II Medical Devices designed by Golden Technologies are intended to assist persons with chronic physical disablements.

Wireless disassembly in 4 easy steps!



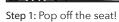



Step 2: Pull out the battery pack!



Step 3: Fold down tiller and lock in place!



Step 4: Pull forward on the drivetrain release lever to separate the front and rear frame sections!



Disassembles into 6 pieces!

4 wheel model shown



A battery charging port is conveniently located on the tiller under the control panel.



Phone: 800-624-6374 Fax: 800-628-5165 www.goldentech.com



^{2.} Electronics warranty excludes batteries.